MAXIWELL

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Manual

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Main voltage: 230V AC
Input: DC 12V
Rated Power: 14.4W

Frequency: 50HZ
Currency: 1.2A
Certified Electrical Standard: CE, CUL, PCT, SAA
Zubehör: 12V adaptor/230V adaptor

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
Thank you for purchasing the Maxiwell. Before you begin to use this massager, please read the operation manual thoroughly for the correct methods of usage, especially the notes on caution. This manual should be kept available for future reference after you have read it through.

Functions and Features

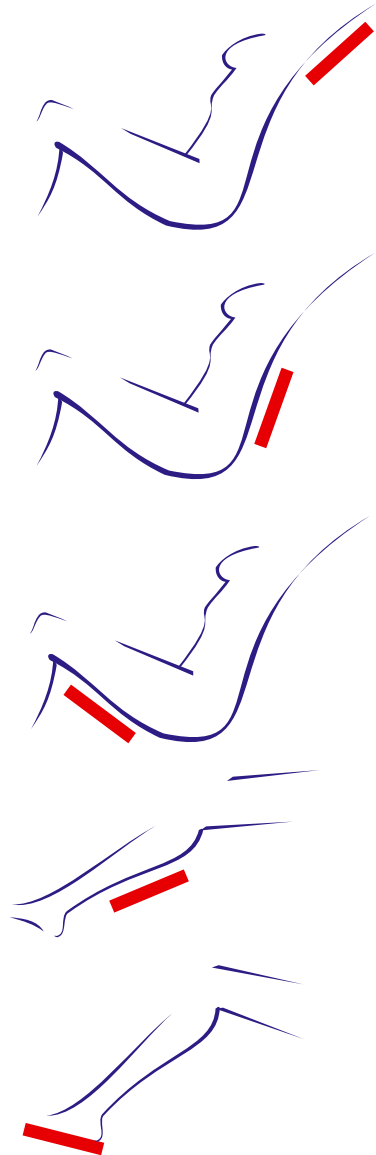4 massage heads massage your back by circular movements

The neck, shoulders and the back can be massaged very comfortably.

Sit down in a chair with a high seat-back and put the Maxiwell like a cushion into the neck.

By slightly bending the upper part of the body forwards the Maxiwell glides during the massage down to the lumbar vertebra. A unique beautiful feeling overcomes you and the strengthening of the back musculature leads in a long term to a correction of the upright body attitude.

Regular and well-performed massage improves blood circulation and -supply, relieves tension in the neck and shoulder region and all together strengthens the muscles of the back. Massaging the back stimulates the nerves to function properly and at the same time the metabolic processes are activated. Massage of the reflex points of the soles of the feet elicits a positive and harmonising effect on almost all of the body organs.

Application duration: The duration of treatment should not exceed 30 minutes. Please avoid directly massages on the spinal column. After a treatment you should wait approx. 15 minutes, before you use the Maxiwell again. Thus you contribute to extend the life span of your Maxiwell.

Reference to the use: it is not recommended to exceed a massage duration of 3 minutes per body area. In case of a too high massage pressure, you can put a towel between the body and the cushion.

New functions

The new Maxiwell massages in both directions – clockwise and counterclockwise. In addition it has 2 speed levels, now – low and high.

Medical references

If your physician/orthopaedist did not expressly forbade the back massage, you can enjoy the massage without doubts.

To be massaged should never be felt as unpleasant. Exaggerated massage can also damage the body.

Pregnant women, patients and children should not use the cushion. With injuries of the spinal column or after operations, a massage can be dangerous. In case of doubts please consult your physician.

Safety references

The Maxiwell is safety examined CE and corresponds to the recognized rules of the technology as well as the valid safety regulations. Despite this passed safety inspections you have to take following points into account:

- Separate the Maxiwell after use from the electric circuit
- Do not pull on the cable of the power supply to unplug it from the socket
- Attach the Maxiwell only to alternating current (see identification plate on the back)
- For maintenance of the Maxiwell please use a soft cloth and a neutral cleaner
- Do not use a solvent containing cleaner
- In case of disturbances, please do not repair yourself. Return the Maxiwell to the responsible customer service agency. With self-repairs any warranty claim expires; be aware of danger.
- Do not use in damp areas, like the bathroom.
- Please never cover the exhaust-openings, because the Maxiwell could overheat.
- Stop the massage immediately, if you feel any pain
- Do not use the cushion while you are driving

Terms of warranty

The dealer grants a legally prescribed guarantee time on the product. The warranty applies starting from the purchase date. The purchase date is to be proven with the purchase voucher.

The dealer will repair in Germany seized products with determined errors, free of charge and/or will replace it. Wearing parts e.g. material covers, are not concerned. A condition for this is the entry of the defective product with purchase voucher before expiration of the guarantee period.

The warranty claim is void, if after statement of an error, e.g. by outside effect or as consequence of repair or change, which was not made by the manufacturer or an authorized contract dealer, appeared.

The warranty granted by the dealer is limited to repair and/or exchange of the product. In the context of this warranty the manufacturer and/or dealer, transfers no adhesion and is not responsible for damage, which results from ignoring the operating instructions and/or abusive application of the product.

If it turns out with examination of the product by the dealer, that the warranty demand concerns an error not seized by the warranty or that the warranty term ran off, the costs of examination and repair is to be carried from the customer.

Technical data

Main voltage: 230V AC

Input: DC 12V

Rated Power: 14.4W

Frequency: 50HZ

Current: 1.2A

Certified Electrical Standard: CE, CUL, PCT, SAA

Accessories: 12V adaptor/230V adaptor

Massage and Health

It's so easy to remain fit and healthy

Regular and well-performed massage improves circulation and blood supply, relieves tension in the neck and shoulder region and altogether strengthens the muscles of the back. In addition, it promotes good overall health and has a soothing effect if you are under stress. With your personalised massage treatment you can ensure that our spine is supported by a strong corset of muscles, thus giving you more years of health, mobility and fitness.

And it is true also for those who are active in sports: Sport and massage go together. You can avoid the pain of a strained or ruptured muscle by a balancing and relaxing massage. It optimises the blood circulation through muscles and tissue and prevent muscular soreness after exercise. Top athletes have always known that.

The Maxiwell supports effectively your personal health program.